


MAGIC MASSAGE PRO

GETTING STARTED

- 1** SNAP ELECTRODE
WIRE to MASSAGE
PADS

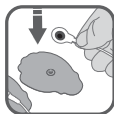

- 2** REMOVE PLASTIC
FILM from PADS



- 3** PLACE PADS ON SKIN
(Refer to the Pad
Placement guide)



- 4** CONNECT ELECTRODE
WIRE to CHANNEL A
*REPEAT 1 - 4 for
CHANNEL B*



5 TURN the unit ON



6 PRESS (M) to CYCLE through THERAPIES



7 PRESS - OR + to ADJUST INTENSITY



8 PRESS the CENTER BUTTON to pause or resume your therapy



9 PRESS the (T) BUTTON to adjust the duration of your therapy

